
# seca 797
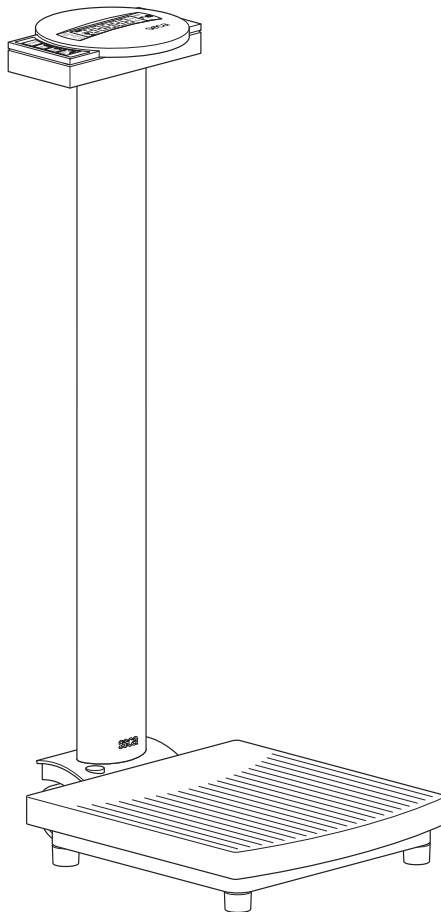
<b>D</b>	<b><i>Bedienungsanleitung und Garantieerklärung</i></b>	<b>3</b>
<b>GB</b>	<b><i>Instruction manual and guarantee</i></b>	<b>15</b>
<b>F</b>	<b><i>Mode d'emploi et garantie</i></b>	<b>27</b>
<b>I</b>	<b><i>Manuale di istruzioni e garanzia</i></b>	<b>39</b>
<b>E</b>	<b><i>Manual de instrucciones y garantía</i></b>	<b>51</b>
<b>DK</b>	<b><i>Betjeningsvejledning og garantibevis</i></b>	<b>63</b>
<b>S</b>	<b><i>Bruksanvisning och garanti</i></b>	<b>75</b>
<b>N</b>	<b><i>Bruksanvisning og garantierklæring</i></b>	<b>87</b>
<b>FIN</b>	<b><i>Käyttöohje ja takuu</i></b>	<b>99</b>
<b>NL</b>	<b><i>Bedieningshandleiding en garantieverklaring</i></b>	<b>111</b>
<b>P</b>	<b><i>Instruções de utilização e declaração de garantia</i></b>	<b>123</b>
<b>GR</b>	<b><i>Οδηγίες χειρισμού και εγγύηση</i></b>	<b>135</b>

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# Full certification

seca products you are purchasing not only technology matured over a century, but also quality confirmed by authorities, the legal system and by institutes.

seca products comply with European directives, standards and national laws. With seca, you are buying the future.

**M**

The products in these operating instructions comply with German legislation governing medical products and therefore with Directive 93/42/EEC of the Council of the European Community which is embodied in national laws throughout Europe.

Scales carrying this symbol are approved for medical purposes in the European Community. seca scales with this symbol meet the stringent quality and technical requirements placed on scales which can be calibrated.

Scales carrying this symbol are calibrated in accordance with Precision Class IIII of EC Directive 90/384/EEC.

Products carrying this symbol comply with the following directives and standards:

1. Directive 90/384/EEC governing non-automatic weighing instruments
2. Directive 93/42/EEC governing medical devices
3. DIN EN 45501 governing metrological aspects of non-automatic scales

seca's professionalism has also been officially recognized. The certificate from TÜV Product Service, the body responsible for medical devices, confirms that as a manufacturer of medical devices, seca rigorously complies with stringent legal requirements. seca's quality assurance system covers the development, production, sales and service of medical scales and measuring systems and ergometry.

seca helps the environment.

We are concerned to save natural resources. We therefore strive to save packaging material where meaningful. What is left over can conveniently be disposed of locally via the German Dual System (recycling programme).

# 1. Congratulations!

---

In the **seca 797** electronic personal scale, you have purchased a high-precision yet sturdy piece of equipment.

seca has been putting its experience at the service of health for over 150 years now, and as market leader in many countries of the world, is always setting new standards with its innovative developments for weighing and measuring.

In accordance with national regulations, the **seca 797** electronic personal scale is used primarily in hospitals, doctors' surgeries and in-patient nursing facilities. The scale is approved to CE class III.

Weight is determined within a few seconds.

In addition to the conventional determination of weight, the **seca 797** also has a function for determining body mass index. For this purpose, height is entered on the keypad and the body mass index associated with the weight value is output.

The scale is of a very sturdy construction and will give you long and loyal service. It is easy to operate and the large display is easy to read.

The **seca 797** can be moved on castors, and its low power consumption means that it can be used on a mobile basis with batteries for long periods.

## 2. Safety

---

Before using the new scale, please take a little time to read the following safety instructions.

- Follow the instructions in the instructions for use.
- Keep the operating instructions and the declaration of conformity in a safe place.
- Do not drop the scale or subject it to violent shocks.
- Use only the type of battery stated (see "Changing the batteries" on page 24).
- Have scale serviced and re-calibrated on a regular basis (see chapter entitled "Maintenance" on page 24).
- Have repairs carried out only by authorised persons.

## 3. Before you start...


---

### Unpacking

- Remove the packaging and place the scale on a firm, level surface.

The scope of supply includes:

- Weighing platform
- Column with display head
- 2 mounting screws
- Compartment with batteries

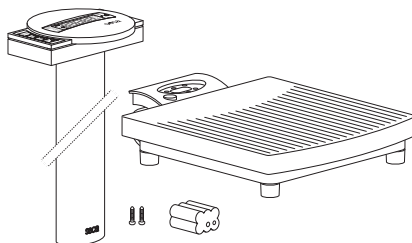
## Power supply

The scale is powered exclusively by batteries. The scale works independently of the mains, which is very handy for mobile use.

The LC display has a low power consumption – 80.000 weighing operations can be performed with one set of batteries. When the battery voltage is no longer adequate, **batt** appears in the display.

- Carefully turn over the scale so that the base of the scale is accessible.
- Push the battery compartment closure towards “Battery” and open the cover.
- Connect the cable to the holder and batteries.
- Place the holder in the battery compartment.
- Close the cover and then turn the scale back up the right way.

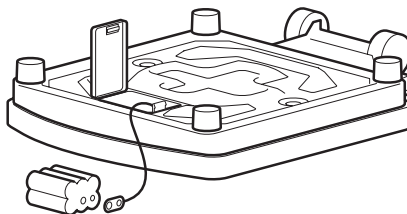
## Fitting the column

The column and display head are fitted to the weighing platform.

- Guide the connecting cable through the opening.
- Place the column in the recess. The display should be legible for people standing on the scale.
- Screw the column to the weighing platform from underneath.
- Connect the plug of the connecting cable to the rear of the scale.
- Then place the scale upright again.

If necessary, lay the scale carefully on its side to tighten up the screws.

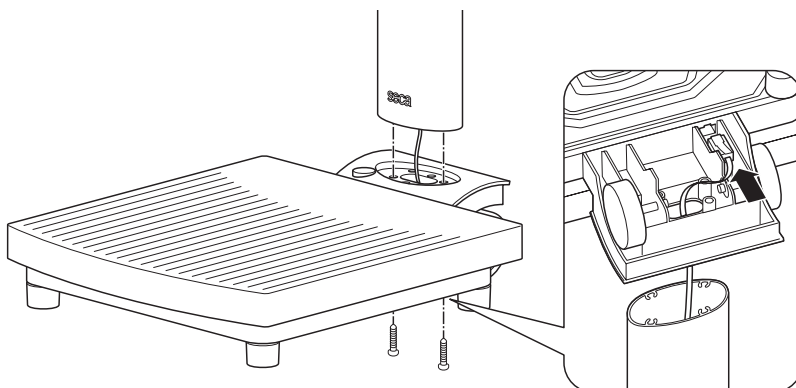
## 4. Setting up and aligning the scale

---

- Place the scale on a level surface.
- Undo the locknuts using the four foot screws.
- Align the scale by undoing the foot screws. The air bubble in the spirit level must be right in the centre of the circle.
- Then re-tighten the locknuts.
- For mobile use on a ward, the **seca 797** can be tilted back slightly and then moved on the castors.

### Important!

**The alignment of the base of the scale must be checked and corrected if necessary every time the location of the scale is changed.**

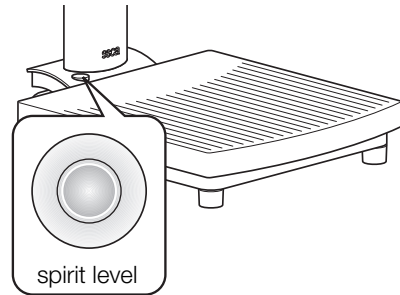
## 5. Operation

---

### Controls and displays

The controls and the display are on the display head on the column.

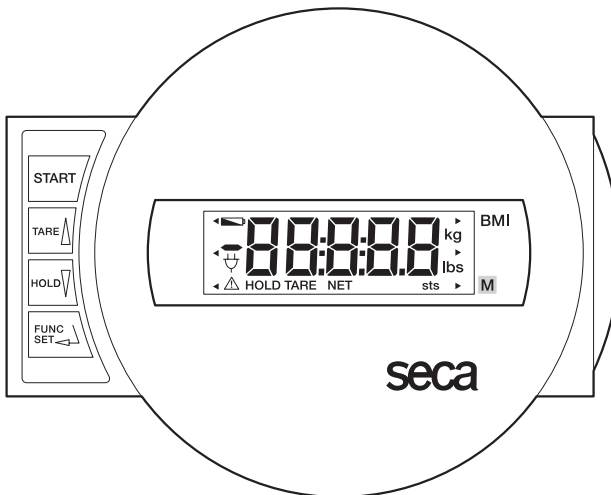Switch scale on and off

press briefly → set display to zero, activate/deactivate tare function  
in function menu → increase value

press briefly → activate hold function  
in function menu → decrease value

Activation of BMI function and Enter key for input data

## Checking that the contents of the calibration counter are correct

This seca scale is calibrated. Calibration may be performed only by authorised centres. To ensure this, the scale is fitted with a calibration counter which records any change to data relevant to calibration. If you want to check whether the scale is properly calibrated, proceed as follows:

- If necessary, switch off the scale.
- Keep any key depressed and start the scale.

The current contents of the calibration counter are shown on the display for a few seconds.

Compare the output contents of the calibration counter with the number shown on the calibration counter sticker.

For valid calibration, both numbers must match. If sticker and calibration counter do not match, then re-calibration is required. Contact your service agent or seca Customer Service.

**Both numbers must match.**

## Correct weighing

- Press the Start key with no load on the scale.  
**SECA**, **88888** and **00** appear consecutively in the display. The scale is then automatically set to zero and ready to use.
- Stand on the scale and keep still.
- Read the weighing result off the digital display.
- If a load greater than 160 kg is placed on the scale, **StOP** appears in the display.
- To switch off the scale, press the Start key again.

## Automatic switch-off/economy mode

The scale switches off automatically after 3 minutes.

## Taring weight display (TARE)

The tare function is very useful, for example for determining the weight of a child, animal or other object placed on the scale in addition to yourself.

Proceed as follows:

- Press the Start key with no load on the scale.

- First stand on the scales yourself without the additional weight and activate the Tare function by pressing the **TARE** key.

The display is reset to zero. **0.0** appears in the display. The note "NET" appears.

- Now take the additional load onto the scale. You can leave the scale to do so, in which case **----** appears.

The scale determines the weight of the additional load.

You can now perform as many weighing operations as you like.

- Press the key down again to deactivate the Tare function; the scale is then back in normal weighing mode and the tare value has been erased.
- The tare value is likewise erased by the scale being switched off.


## Store weight value (HOLD)

You can also continue displaying the weight determined once the load has been removed. This makes it possible to look after a patient first, before having to note the weight.

- Press the **HOLD** key whilst the load is on the scale.

"HOLD" appears in the digital display. The value is stored in the display until the scale switches off automatically after 3 minutes. The triangle with the exclamation mark indicates that this frozen display is an additional facility which is not subject to calibration.

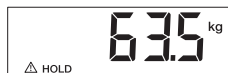- The HOLD function can be switched off by pressing the **HOLD** key again.

00 appears in the display if the scale has no load on it, otherwise the current weight. The "HOLD" display vanishes.

**Note:**

**The HOLD and TARE functions can also be used simultaneously.**

## Body Mass Index (BMI)

Body Mass Index puts height and weight in relationship to one another, providing more accurate information than ideal Broca weight, for example. A tolerance range is given which is considered ideal for health. The Body Mass Index applies only to adults: these values are not relevant to children and young people.

- Press the Start key with no load on the scale.

SECR, 88888 and 00 appear consecutively in the display. The scale is then automatically set to zero and ready for operation.

- Press the **FUNC** key.  
On the display, you will see the last height entered.
- You can change the value in 1 cm increments using the arrow keys.

- Once the correct value is set, press the **FUNC** key again. The BMI function is now activated
- Step on the scale and keep still.
- Read off the BMI on the digital display and compare it with the categories on page 22.
- Press an arrow key if you would like to see the weight briefly.

- To return to normal weighing mode, press the **FUNC** key.
- To switch off the scale, press the Start key again.

## Evaluating Body Mass Index:

Compare the value determined with the categories quoted below which correspond to those used by WHO, 2000.

- **BMI below 18.5:**  
Warning, this patient weighs too little. There could be a tendency to anorexia. An increase in weight is recommended to improve well-being and performance. If in doubt, consult a specialist.
- **BMI between 18.5 and 24.9:**  
The patient is of a normal weight. He can stay as he is.
- **BMI between 25 and 30 (pre-obese):**  
The patient is slightly to moderately overweight. He should cut his weight if he is already ill (e.g. diabetes, hypertension, gout, fat metabolism disorders).
- **BMI over 30:**  
Weight reduction is essential. Metabolism, cardiovascular system and bones are all under strain. A consistent diet, plenty of movement and behavioural training are all recommended. If in doubt, consult a specialist.

## 6. Cleaning

---

Clean the rubber surface and the housing as required using a domestic cleaning agent or commercially-available disinfectant. Follow the manufacturer's instructions.

Under no circumstances use abrasive or acid cleaners, white spirit, benzene or the like for cleaning. Such substances can damage the high-quality surfaces.

## 7. What do I do if...

---

### ... no weight display comes on?

- Is the scale switched on?
- Check the batteries.
- Is the plug contact on the scale properly connected??

### ... 0.0 does not appear before the weighing operation?

- Press the Start key again – there must not be any load on the scale

### ... - - - - appears instead of 0.0 before weighing?

- Press the Start key again – there must not be any load on the scale.

### ... one segment is illuminated constantly or not at all?

- The relevant segment has a fault. Inform Servicing.

### ... the display shows ?

- Battery voltage is running low. Change the batteries in the next few days.

### ... bAtE appears in the display?

- Batteries are empty. Put in new batteries.

### ... 5t0P appears in the display?

- Maximum load has been exceeded.

### ... the display flashes?

- If you have just activated the Hold or Tare function, please wait until the scale has detected a steady weight value. The display will then stop flashing.
- If there has been no activation of a function beforehand, remove the load from the scale and wait until 0.0 is displayed, then weigh again.

### ... the display tENP appears?

- The ambient temperature of the scale is too high or too low. Place the scale in an ambient temperature between +10 °C and +40 °C. Wait about 15 minutes for the scale to adapt to the ambient temperature and then weigh again.

### ... the display E and a number appear?

- Switch off the scale with the Start key and start the scale again. The scale will then work normally again.

If this is not the case, disconnect the power supply by briefly removing the batteries. If this measure is equally unsuccessful, inform the Service department.

## 8. Maintenance

---


Have the scale recalibrated by authorised persons in accordance with national legal requirements. The year of initial calibration can be found behind the CE symbol on the rating plate above the number of the appointed office, 0109 (Hessen Calibration Office). Recalibration is essential if one or more security marks are damaged or the calibration counter does not match the number on the applicable calibration counter mark (cf. page 19).

We recommend having your scale serviced by your local service agent before having it recalibrated. seca Customer Service will be pleased to give you further assistance.

### Changing the batteries

When  or **bAtE** appear in the display, you should change the batteries. You need four type AA 1.5 V batteries.

- Carefully turn over the scale so that the base of the scale is accessible.
- Push the battery compartment closure towards “Battery” and open the cover.
- Take out the holder and batteries.
- Remove the used batteries and insert new batteries in the holder.
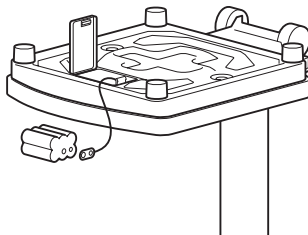
### Check that the polarity is correct.

- Place the holder and batteries in the battery compartment.
- Close the cover.
- Then turn the scale back up the right way.

## 9. Technical data

---

### Dimensions

width:	292 mm
height:	830 <sup>+10</sup> mm
depth:	390 mm

### Height of figures

20 mm

### Weight

approx. 6.8 kg

### Temperature range

+10 °C to +40 °C

### Batteries

4 type AA 1.5 V batteries

### Calibration as per Directive 90/384/EEC

approved for medical use, Class III

### Medical product as per Directive 93/42/EEC

Class I

### Colour

white

### Maximum load

160 kg

### Minimum load

2 kg

### Graduations

200 g

### Accuracy at first calibration

0 to 10 kg:	± 100 g
10 kg to 40 kg:	± 200 g
40 kg to 160 kg:	± 300 g

## 10. Spare parts and accessories

---

### Measuring rods:

material: hard-anodised aluminium, slide with lockable measuring tab

### seca 220 telescopic measuring rod

seca model 220-17-14-004:

Measuring range: 820 - 2000 mm,

Graduations: 1 mm

seca model 220-18-14-004

(double graduations):

Measuring range: 820 - 2000 mm and

32 1/4 " - 78 3/4 ",

Graduations: 1 mm and 1/8 "

# 11. Disposal

---

## Batteries

Do not throw used batteries away in domestic waste. Dispose of batteries at collection points in the vicinity. When buying new batteries, select those low in harmful substances and containing no mercury (Hg), cadmium (Cd) or lead (Pb).

## Disposing of the scale

If the scale can no longer be used, the responsible waste disposal association will be glad to advise you on the measures necessary to dispose of it properly.

# 12. Warranty

---

A two-year warranty from date of delivery applies to defects attributable to poor materials or workmanship. All moveable parts - batteries, cables, mains units, rechargeable batteries etc. - are excluded. Defects which come under warranty will be made good for the customer at no charge on production of the receipt. No further claims can be entertained. The costs of transport in both directions will be borne by the customer should the equipment be located anywhere other than the customer's premises. In the event of transport damage, claims under warranty can be

honoured only if the complete original packaging was used for any transport and the scale secured and attached in that packaging just as it was when originally packed. All the packaging should therefore be retained.

A claim under warranty will not be honoured if the equipment is opened by persons not expressly authorised by seca to do so.

We would ask our customers abroad to contact their local sales agent in the event of a warranty matter.

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